

# STARFORMER™

Magnetic Body Sculpting and Muscle Strengthening



Fotona  
choose perfection

# STARFORMER™

StarFormer™ is Fotona's latest innovation in body shaping. This breakthrough technology utilizes magnetic pulse energy to stimulate muscle contraction. As a result, muscles strengthen, tone and firm.


The system provides two treatment options:

- **TightWave™**  
muscle strengthening, toning and firming of weakened areas
- **IntimaWave™**  
strengthening of pelvic floor muscles to treat stress urinary incontinence



## How it Works

StarFormer™ uses magnetic stimulation technology to selectively target muscle tissues. Magnetic pulses harmlessly penetrate into the muscle layers causing them to automatically contract. The end result is increased muscle mass.

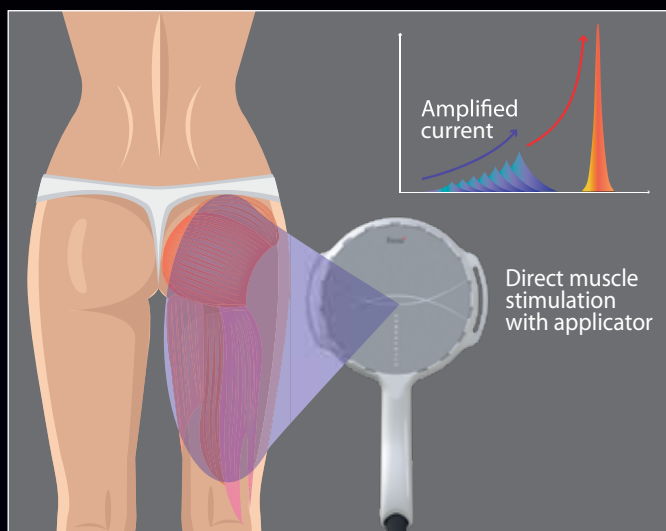

# TIGHTWAVE™

TightWave™ is a breakthrough treatment for effectively strengthening, shaping and firming different areas of the body:



- Glutes
- Abdomen
- Hamstring area
- Arms
- Back muscles




The treatment is performed with the patient lying down fully clothed. One session generally lasts between 20 to 30 minutes with patients resuming their daily activities immediately after treatment.

# INTIMAWAVE™

StarFormer™ can additionally be used to stimulate pelvic floor muscles and back, utilizing two applicators to treat all types of male or female urinary and faecal incontinence. During treatment the patient remains fully clothed and sits on a comfortable IntimaWave™ chair. The magnetic applicators in the chair stimulate muscle contraction in the entire pelvic floor area, increasing the strength and endurance of the pelvic floor. The result is effective and long-term relief of uncontrolled urine flow.




Before

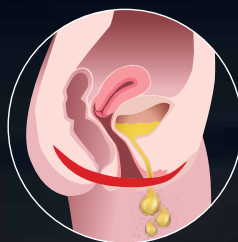

*Pelvic floor muscles insufficiently support pelvic organs effecting bladder control.*

After

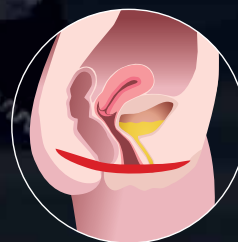


*Magnetic stimulation induces pelvic floor muscle strengthening and the return of normal continence function.*



# One Device, Multiple Solutions

## Key Benefits

- Safe, quick and easy muscle strengthening
- Non-invasive
- Improvement after only a few sessions
- Walk-in/walk-out procedure
- High success rate and patient satisfaction

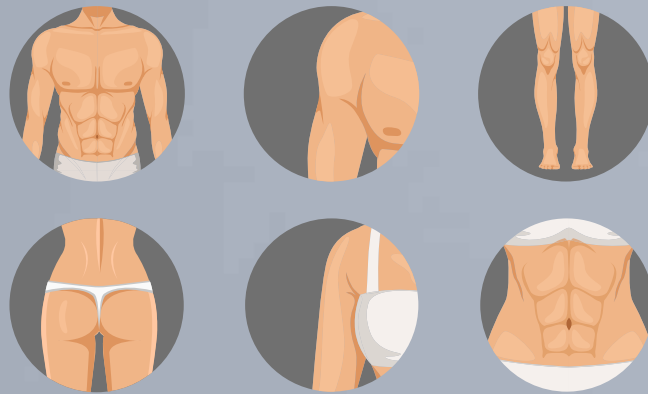

|                               |                          |
|-------------------------------|--------------------------|
| Device name:                  | Fotona StarFormer™       |
| Input voltage:                | 100-240V AC, 50~60Hz     |
| Power:                        | 1500 VA                  |
| IP protection:                | IPX0                     |
| Size (W x L x H):             | 421 mm x 630 mm x 843 mm |
| Weight:                       | 50 kg                    |
| IEC 60601/1 classification:   | Class I, Type BF         |
| MDD 93/42/EEC classification: | Class IIa                |

**CE**  
2265

All Fotona medical devices are CE marked and approved to be sold in the EU. For countries where specific national approvals or clearances are required, some of the products and/or applications may not yet have been approved. Please check with Fotona, your local Fotona distributor or your national regulatory body about whether a specific product or application has been approved to be marketed and sold in your country.

For related patents see: [www.fotona.com/patents](http://www.fotona.com/patents)

[www.fotona.com](http://www.fotona.com)

Fotona, d. o. o.  
Stegne 7  
1000 Ljubljana  
Slovenia  
EU

Fotona, LLC  
2307 Springlake Road  
#518  
Dallas, TX 75234  
USA

Fotona GmbH  
Hohlbachweg 2  
73344 Gruibingen,  
Germany  
EU

Fotona Beauty Light, (Suzhou)  
Medical Devices Co, Ltd.  
No 2, Zengfu Road, Guli Town  
Changshu City, Jiangsu Province  
CHINA, 215515